


## PROVING RINGS

The **G3831** thru **G3844** are single compression proving rings with dial gauge reading to .0001". The rings are constructed from heat treated, plated, alloy steel. The blocks are plated steel with a 3/4" tapped hole. Also furnished with the proving rings is a dial gauge and calibration chart. The **G3865** thru **G3872** are double compression proving rings. These proving rings have a high sensitivity at low ends since only one ring deforms until contacting the second ring. The dial gauge has a reading to .001". The rings are constructed from heat treated, plated, alloy steel. The blocks are aluminum with 1/4" tapped hole. Also furnished with dial gauge and calibration chart.

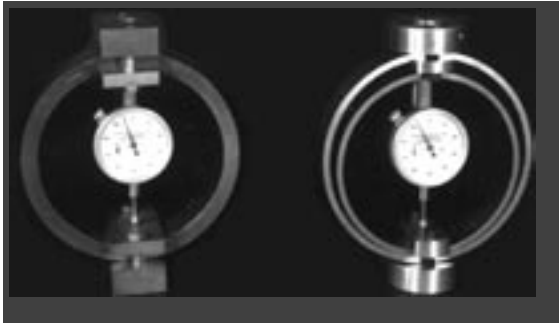

MODELS	DESCRIPTION	LOAD CAPACITY
<b>G3831</b>	SINGLE PROVING RING	1,000 LBS
<b>G3832</b>		2,000 LBS
<b>G3836</b>		6,000 LBS
<b>G3841</b>		10,000 LBS
<b>G3843</b>		15,000 LBS
<b>G3844</b>		20,000 LBS
<b>G3863</b>	DOUBLE PROVING RING	300 LBS
<b>G3865</b>		500 LBS
<b>G3867</b>		750 LBS
<b>G3871</b>		1,000 LBS
<b>G3872</b>		1,500 LBS

## OHAUS MECHANICAL BALANCES



These balances combine a history of outstanding engineering and design expertise which have made them classic weighing instruments. Some of the outstanding features of these instruments are notched and tiered beams, positive poise positioning, specific gravity weighing, and magnetic dampening. Some of the other features include covered agate bearings, out of level accuracy, and precision ground knives. Specific gravity weighing is possible on some models without any accessories, however, some of the balances need an extra rod and clamp assembly (**G4189** sold separately).

MODELS	DESCRIPTION	CAPACITY	READABILITY	PLATFORM	SPECIFIC GRAVITY
<b>G4121</b>	CENT-O-GRAM BALANCE	311 G	0.01 G	3.5 X 0.5	YES
<b>G4131</b>	DIAL-O-GRAM BALANCE	310 G	0.01 G		YES
<b>G4170</b>		610 G / 2610 G	0.1 G / 0.01 G	6" DIA.	W/ <b>G4189</b>
<b>G4171</b>		2610 G	0.01 G		W/ <b>G4189</b>

ACCESSORIES	DESCRIPTION
<b>G4134</b>	DUST COVER
<b>G4135</b>	CARRYING CASE FOR <b>G4121</b> & <b>G4131</b>
<b>G4174</b>	FOOTED SCOOP, ST. STL., 12" X 6" X 2.75"
<b>G4185</b>	EXTRA WEIGHT, 500 G
<b>G4186</b>	EXTRA WEIGHT, 1,000 G
<b>G4187</b>	SCOOP, ST. STL., 12" X 6" X 2.75" W/ WEIGHT
<b>G4189</b>	CLAMP & ROD ATTACHMENT SET
<b>G4190</b>	SUPER SCOOP